

**PAJUNK®**

## *Vario Balloon Systems*

*The safe solution with universal seal*




Balloon systems

# PAJUNK® Vario Balloon Systems NEW

## The safe solution with universal seal

Universally applicable access and anchoring systems are growing in importance in open laparoscopy, hernia and bariatric surgery. The new PAJUNK® Vario balloon systems with working channel are equipped with a universal double seal system. They are available as structural, ring-anchor and URO ring-anchor balloon systems.

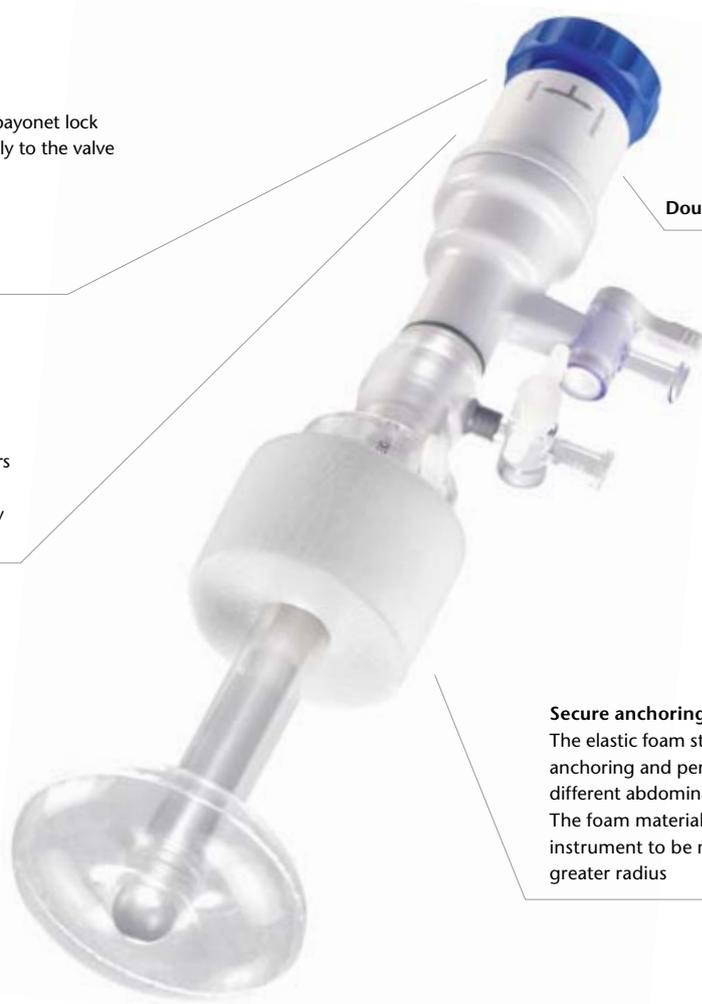



**Lockable obturator**  
When inserted completely, a bayonet lock connects the obturator securely to the valve housing.

**Simple unlocking**  
A twist is all that's needed.



**Universal seal**  
All instruments with diameters between 5 and 11 mm can be inserted with a completely gas-tight seal.



**Secure anchoring**  
The elastic foam stopper offers secure anchoring and permits adaptation to different abdominal wall thicknesses. The foam material also allows the instrument to be maneuvered in a greater radius.

**Three balloon variants**



Structural balloon



Ring-anchor balloon



URO ring-anchor balloon

### Vario balloon systems with universal seal

#### Balloons for securing the working access

	Product	Size	Art. No.	PU	Ball pump	Syringe 30ml
	Structural balloon system	Ø 11 mm	1283-20-11S	5	•	
	Ring-anchor balloon system	Ø 11 mm	1283-30-11S	5		•
	Ring-anchor balloon system, short	Ø 11 mm	1283-40-11S	5		•
	URO ring-anchor balloon system	Ø 11 mm	1283-70-11S	5		•
	URO ring-anchor balloon system, short	Ø 11 mm	1283-75-11S	5		•

#### PAJUNK GmbH

Medizintechnologie

Karl-Hall-Strasse 1

D-78187 Geisingen / Germany

Phone +49 (0) 77 04/92 91-0

Fax +49 (0) 77 04/92 91-6 00

[www.pajunk.com](http://www.pajunk.com)